


### 30.10.16a APPENDIX A: Intradialytic Parenteral Nutrition

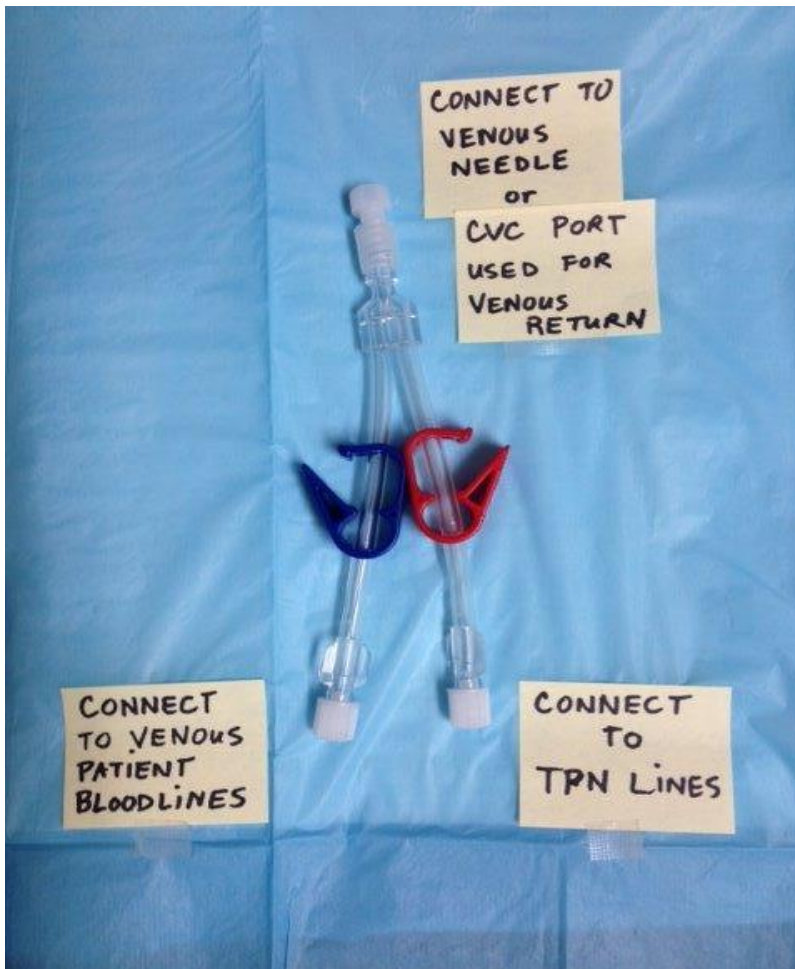

- Prime each arm of the Y-connector with 0.9% Normal Saline pre-filled syringe prior to attaching it to the Central Line venous access.
- If using an AVF, may prime with 0.9% Normal Saline prior to attaching to the venous needle or alternatively, may attach to the venous needle and then backfill with blood by loosening the caps slightly to allow blood to fill. Tighten the caps once filled and flush with 0.9% Normal Saline pre-filled syringe.
- When ready to initiate Hemodialysis, attach venous bloodline to blue lumen of y-connector and attach TPN lines to red lumen of y-connector.
- Attach arterial bloodline to arterial needle or CVC port used for supply.